

Adaptaid⁺

Transport Cart for Probes and Endoscopes



Endoscopy
Sterilization/MDR

Ultrasound
Operating Room

Probe and Endoscope Transport Carts

Essential for the safe transport of probes and flexible endoscopes, these carts have a robust structure made entirely of grade 304 stainless steel. Thanks to an ergonomic wide grip handle and stainless steel wheels with thermoplastic rubber treads, movement is effortless and silent.

The cart interior has sliding rails with retaining stops to keep the trays in place when moving. Five to six trays with lids are included, depending on the cart model selected. The carts and trays are compatible with the automated washers used in hospitals.

MODELS	DIMENSIONS			STRUCTURE	DOORS	TRAYS
	INCHES	W X D X H	MM			
EDS-030-P005-O	28.5 x 20.3 x 38.5	724 x 514 x 978		open on 2 sides	0	5
EDS-030-P005	28.5 x 20.3 x 38.5	724 x 514 x 978		open on 1 side	0	5
EDS-030-P005-F	28.5 x 20.3 x 38.5	724 x 514 x 978		closed	2	5
EDS-030-P006-O	28.5 x 20.3 x 43.5	724 x 514 x 1105		open on 2 sides	0	6
EDS-030-P006	28.5 x 20.3 x 43.5	724 x 514 x 1105		open on 1 side	0	6
EDS-030-P006-F	28.5 x 20.3 x 43.5	724 x 514 x 1105		closed	2	6



Light

EDS-030-P006-O

Ideal for tight spaces, this cart provides quick access to trays from both the front and rear.



Practical

EDS-030-P006

The rear panel and front opening provide a user-friendly solution that combines tray protection and accessibility.



Efficient

EDS-030-P006-F


Equipped with two doors, this cart offers increased protection for the trays it transports.

Transport Cart Features

Quality Components for Durability, Ergonomics and Productivity

Structure	Solid structure with 304 #4 stainless steel assembly. Lightweight frame, easy to clean and compatible with cart washers.
Slide	Made of stainless steel with retaining stops for safe movement.
Doors*	Gasket sealed. Drainage holes. Locking device to secure doors in the open position (270° opening).
Handle	Made of stainless steel with a diameter of 1.5 inches (40 mm), as recommended by ergonomists. Can be installed on the left or right side (as required).
Casters	Five-inch (125 mm) diameter for the four swivel and multidirectional stainless steel casters. Precision ball bearing, wire guard and thermoplastic rubber band. Total brake on two casters.
Accessories	Trays with transparent lids and reversible, reusable “clean/dirty” labels.

* Available only with model EDS-030-P001-F.

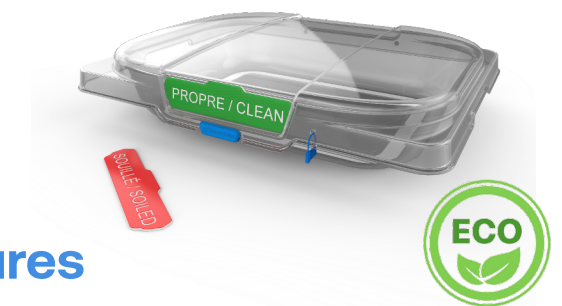

Probe and Scope Tray Features

Compatible with Large Probes and Flexible Endoscopes

Structure	Molded case in transparent gray polycarbonate with two side handles and a built-in hole to insert tamper-evident locks. Robust, durable and environmentally friendly. No consumables required.
Lid	Molded in transparent polycarbonate designed for quick viewing of content.
Dimensions	22.9 W x 18 D x 5 H inches, 582 W x 457 D x 127 H mm
Cleaning	Easy to clean and compatible with automated washers.

Adaptaid⁺

High-Performance Medical Equipment
Designed for Healthcare Professionals

For more information about our products, contact our customer service department.



CONTACT US
+1 450-585-7272
customerservice@adaptaid.com



©2026 ADAPTAIDE INC. ALL RIGHTS RESERVED

adaptaid.com