

Adaptaid<sup>+</sup>

---

# ERGOFLEX

## Prep & Pack Workstation



Sterilization/MDR

## ERGOFLEX Workstation

Developed to meet the diverse needs of healthcare workers, the ERGOFLEX workstation provides an ergonomic, comfortable workspace for easy handling of instruments to be reprocessed.

In addition to ergonomics and safety for healthcare workers, its design also focuses on productivity, offering a wide range of sizes and accessories to ensure a configuration that meets specific user requirements in every respect.

## A Multifunctional and Efficient Workstation

ERGOFLEX workstations are designed to facilitate the inspection, assembly, and packaging of surgical instruments.

To create a customized configuration, the table has around twenty accessories that can be added to the worktop or to the upper module, which increases storage potential and provides optimal lighting.

MODELS	WORKTOP DIMENSIONS			DIMENSIONS WITH UPPER MODULE		
	INCHES	W X D	MM	INCHES	W X D X H	MM
RDM-010-P002-EX6030	60 x 30		1524 x 762	60 x 34 x 76		1524 x 864 x 1930
RDM-010-P002-EX6036	60 x 36		1524 x 914	60 x 40 x 76		1524 x 1016 x 1930
RDM-010-P002-EX7230	72 x 30		1829 x 762	72 x 34 x 76		1829 x 864 x 1930
RDM-010-P002-EX7236	72 x 36		1829 x 914	72 x 40 x 76		1829 x 1016 x 1930
RDM-010-P002-EX8436	84 x 36		2134 x 914	84 x 40 x 76		2134 x 1016 x 1930



### Compact

RDM-010-P002-EX6030  
Ideal for small spaces, this table offers plenty of storage space on the lower shelf.




### Functional

RDM-010-P002-EX7230  
A user-friendly solution that combines ergonomics and accessibility thanks to the upper module and proximity of accessories.



### Efficient

RDM-010-P002-EX7236 & RDM-010-SC  
Equipped with a sterile wrap rack and an upper module, this package creates a complete, ergonomic workstation.

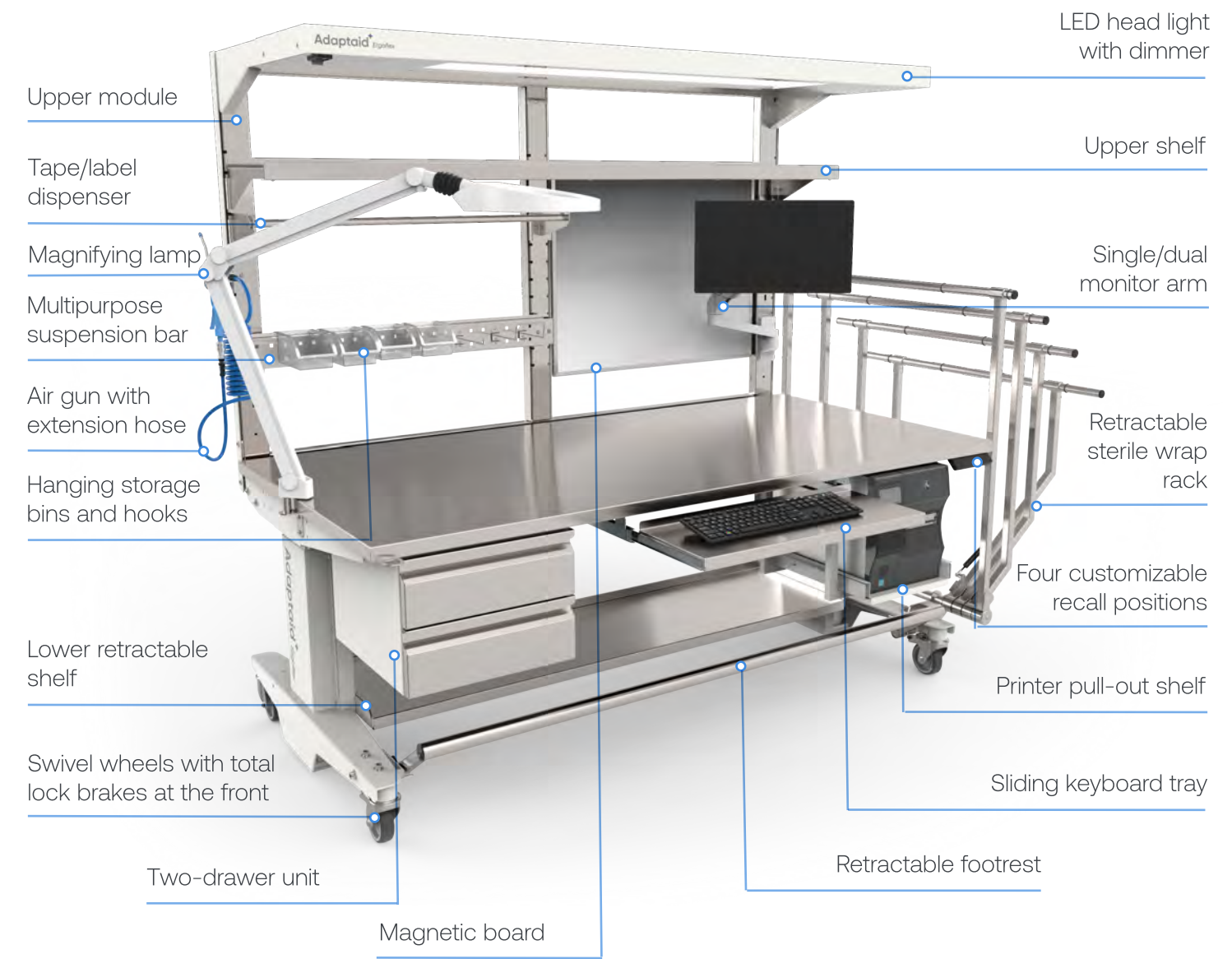

## Retractable Sterile Wrap Rack

Installed directly on the base of the ERGOFLEX preparation and packing table, this retractable sterile wrap rack optimizes the storage of sterile wraps by reducing the space required on work surfaces, walls and floors.



## An Ergonomic, User-Friendly, and Productive Solution

The ERGOFLEX retractable sterile wrap rack provides fast, close and ergonomic access to sterile wraps.

What's more, its proximity to the table increases the comfort and efficiency of the sterile processing staff in carrying out the packaging process of medical devices. Among its many benefits are:

- + Reduced footprint thanks to the retractable frame.
- + Optimization of available floor space.
- + Ergonomic organization of sterile wraps.
- + Reduced musculoskeletal injuries thanks to lateral access to sterile wraps.
- + Increase productivity by enabling a continuous workflow.
- + Compatible with various sizes of sterile wraps.


### Solid Structure

Robust, torsion-resistant 304 stainless steel structure. Four cylindrical bars hold the sterile wraps in place with retaining rings.



#### OPEN STERILE WRAP RACK

- + Quick and effortless opening
- + Direct access to wraps
- + Depth of 19 7/8 inches (505 mm) when open

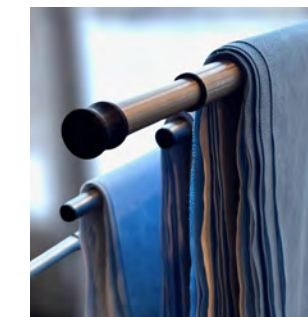

### Additional Options

- + Telescopic extension bar 40.4 to 51.5 inches (1026 to 1308 mm)
- + Fixed, non-retractable frame (model RDM-010-SCF)

#### CLOSED STERILE WRAP RACK

- + Closes with a single push of the hand
- + Increased lateral clearance
- + Depth of 12 inches (305 mm) when closed



MODELS	DIMENSIONS			*SPACE SAVING	
	INCHES	W X D X H	MM	INCHES	MM
RDM-010-SC	39.6 x 19.8 x 45.4	1006 x 505 x 1228		7.8	198
RDM-010-SCF	39.6 x 20 x 45.4	1006 x 508 x 1228		N/A	N/A

\* Lateral space gained when the retractable rack is closed.



## High-Performance Standard Features

Designed with comfort and safety in mind, ERGOFLEX workstations are ideal for a whole range of tasks relating to the reprocessing of medical devices.



### Mobile and Robust Table

The baked-on, electrostatic-coated steel table bases have been carefully designed to provide a sturdy base perfectly suited for integrating the optional ERGOFLEX retractable sterile wrap rack. Boasting excellent strength, the table's load capacity is 750 lb (340 kg).

For ease of movement, the legs are equipped with four swivel casters. The two front casters are equipped with a total brake that can be activated in the stationary position with a simple press of the foot.



### The Ideal Worktop

Resistant to hospital disinfectants and repetitive cleaning, the non-porous worktop is made of high-quality, low reflection grade 304 stainless steel with an antibacterial sealed underside.

## Configurable Accessories to Enhance Ergonomics, Productivity and Safety

ERGOFLEX workstations are built for ergonomics, productivity and user safety, with a wide range of accessories to optimize instrument inspection, preparation and packaging.

### Quick Height Adjustment

Each model features an electrically height-adjustable worktop. The height of the worktop can vary from 30 to 45 inches (762 to 1143 mm), depending on the position required by the user. It is also possible to store four different positions in memory to provide quick height adjustment for multiple users of the same workstation.



The sliding tray for the keyboard and mouse can be put back in its place as soon as it is no longer needed.

- + CPU holder
- + Drawer unit
- + Magnifying lamp with adjustable arm
- + Printer holder
- + Retractable footrest
- + Sliding keyboard tray
- + Sliding lower shelf
- + Sterile wrap rack
- + Tabletop under-mount for power strip

The height adjustment makes it easy to switch between standing and sitting positions.



## Useful Accessories for the Upper Module

Perfectly integrated with the ERGOFLEX table, the top module follows the up and down movements of the worktop to provide the optimum height for each user. The top module also acts as a storage space, keeping essential reprocessing supplies and accessories within easy reach.

- + Air gun with holder
- + Inclined storage shelf
- + LED lighting head with dimmer
- + Magnetic whiteboard
- + Multipurpose suspension bar
- + Single/dual monitor arm
- + Storage bins (various sizes)
- + Tape/label dispenser
- + Upper storage shelf







## Storage Bins

This polycarbonat bin is autoclavable and comes with a transparent lid for quick and easy identification of contents. The lid can be locked in the open position for easy access to supplies.

## Autoclavable Containers, Useful and Durable

- + The sturdy molded edge is easily attached to the suspension bar of the ERGOFLEX workstation's upper module.
- + Four sizes available with optional lids for greater storage versatility.



OVERVIEW	BIN MODELS	OVERALL DIMENSIONS			LID MODELS
		INCHES	W X D X H	MM	
	ADA-01610	5 3/8 x 4 1/8 x 3 1/4	137 x 105 x 83	ADA-01611	
	ADA-01612	7 3/8 x 4 1/8 x 3 1/4	187 x 105 x 83	ADA-01613	
	ADA-01614	10 7/8 x 4 1/8 x 3 1/4	276 x 105 x 83	ADA-01615	
	ADA-01616	10 7/8 x 5 1/2 x 5 1/4	276 x 140 x 133	ADA-01617	

## Ergonomic accessories

Each ERGOFLEX table has a non-reflective 304-grade stainless steel worktop with an antibacterial sealed underside. The base, with its oven-baked electrostatic-coated steel finish, allows for 15-inch (381 mm) height adjustment via electric control with four customizable memory positions. The load capacity is 750 lb (340 kg) evenly distributed. The base has four 4-inch (102 mm) swivel casters, two of which are equipped with total brakes.



Retractable (or fixed) stainless steel sterile wrap rack



Baked-on, electrostatic-coated steel drawer unit with optional lock (choice of 1, 2, or 3 stainless steel drawer)



Single/dual monitor arm



Stainless steel sliding lower shelf



Tabletop under-mount for power strip with a choice of 8 or 12 outlets



Baked-on, electrostatic-coated steel printer retractable shelf



Magnifying lamp with adjustable arm



Stainless steel sliding keyboard tray



Stainless steel side shelf



Full-length stainless steel footrest



## Accessories for the Upper Module

The ERGOFLEX upper module stands out for its sturdy steel structure with an electrostatic baked-on finish. Available in sizes of 60, 72, and 84 inches (1524, 1829, and 2134 mm), its width corresponds exactly to that of the table on which the upper module is installed.

The stainless steel vertical posts allow accessories to be installed at the desired height in two-inch (51 mm) increments.



LED lighting head with dimmer



Storage shelf with stainless steel anti-slip edge



Stainless steel inclined storage shelf



Stainless steel tape/label dispenser



Adjustable monitor arm with keyboard tray



Articulated arm for single or dual monitors




Magnetic whiteboard



Air gun with hose and holder



Multipurpose stainless steel suspension bar

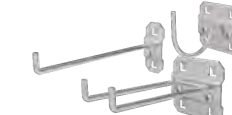

Single, double, and curved hooks in various sizes made of stainless steel

\* The availability of certain additional options may vary depending on the table selected.

# Adaptaid<sup>+</sup>

High-Performance Medical Equipment  
Designed for Healthcare Professionals

**For more information about our products, contact our customer service department.**



CONTACT US  
+1 450-585-7272  
customerservice@adaptaid.com



©2026 ADAPTAIDE INC. ALL RIGHTS RESERVED

[adaptaid.com](https://adaptaid.com)